

25L DIGITAL HEATABLE WATERBATH

Model: 416



FEATURES

- Ideally suited for water or oil applications
- For applications requiring a high degree of accuracy
- Inner tank is manufactured from grade 304 stainless steel for exceptional corrosion resistance
- Reliable, high quality incoloy heating elements
- Heaters and sensors protected by a perforated removable diffuser to ensure even temperature distribution
- Auxiliary heating element to assist in rapid heat up of bath if required
- High density mineral wool insulation for minimum heat loss
- Unit designed for operation at ambient +10°C to 199°C
- Electronic P.I.D temperature control with digital indication of set and process value
- High accuracy pt 100 temperature sensor

NOTE

- Grade 304 stainless steel optional lid



ALL SCIENTIFIC EQUIPMENT COMES WITH A TWO
YEAR WARRANTY ON COMPONENTS AND
DEFECTIVE WORKMANSHIP - Ts & Cs APPLY.

www.scientific.co.za

25L DIGITAL HEATABLE WATERBATH




Model: 416

SPECIFICATIONS

POWER	2000W
VOLTAGE	230V / 50Hz
PACKAGING (W x D x H)	925 x 425 x 330 (mm)
UNIT (W x D x H)	790 x 350 x 250 (mm)
NET WEIGHT	15kg
GROSS WEIGHT	16kg
SHIPPING VOLUME	0.13m ³
TANK DIMENSIONS	460 x 305 x 200 (mm)
TANK VOLUME	25ℓ

OPTIONAL EXTRAS

CODE	DESCRIPTION
133	Stainless steel Gabled lid

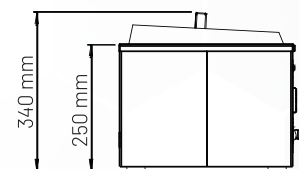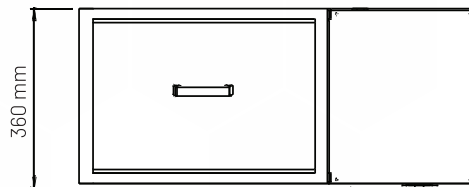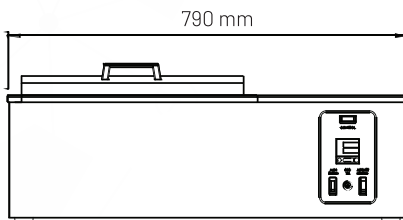

NOTE

- When viewing the appliance from the front in its usual operating position, the width (W) of the product is the total distance from left to right; the depth (D) of the product is the total distance from the front to the back; the height (H) is the total distance from the bottom of the product to the top. Pictures may vary from actual product.

DISCLAIMER

- Output performance figures quoted are dependant on various factors. Scientific reserves the right, without notice, to make changes and revisions to product specifications, materials and design, when we believe it will provide better performance, durability, and efficiency. Unit measurements may vary depending on point of reference.

WARRANTY

- All Scientific equipment comes with a standard ONE year warranty, The 1 Year Extended Warranty is only valid for customers who have registered their warranty online at www.scientific.co.za within the first year after making their purchase.